|EN

AL+ COSMETICS

LumiTech+ Advanced

Internal Observational Skin Follow-Up Report

Internal skin metric follow-up using Aura 3D Face Analysis


Baseline: 13.10.2025	Follow-up: 20.02.2026	Period: approx. 4 months
----------------------	-----------------------	--------------------------

Document purpose

This document presents an internal observational follow-up of selected skin metrics after a period of LumiTech+ Advanced use as part of a regular home skin routine. The data are extracted from two Aura Face Analysis measurements performed in the same clinical setting.

Important: This is not a randomized clinical trial and should not be interpreted as a guarantee of individual results.

NEUTRAL

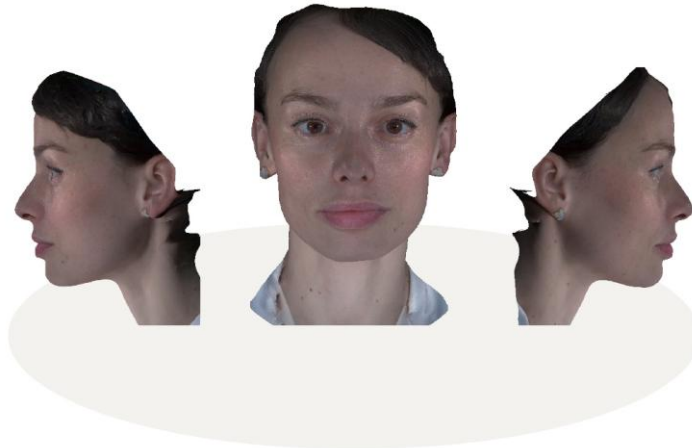Capture Date: 20.02.2026 13:07

Visual excerpt from Aura Face Analysis - facial positioning and analysis views. Personal identifiers from the original reports are not reproduced in this public document.

Prepared for transparency and product education. Version 1.0 / May 2026

Executive summary

The follow-up compares a baseline measurement from 13.10.2025 with a follow-up measurement from 20.02.2026. The Aura system uses a 1-5 scale, where a lower score is more favorable, as well as a relative +/- score compared with people of the same age, gender and skin type.

1.6 -> 1.3 Texture lower score	2.7 -> 2.3 Brown spots lower score	1.2 -> 1.1 Wrinkles lower score	0.0 -> -0.2 Skin score relative score
---	---	--	--

<p>What improved</p> <p>The clearest directional improvements were observed in texture and brown spots. Wrinkles showed a small positive change, while red areas and pores remained stable.</p>	<p>What this means</p> <p>The results support the positioning of LumiTech+ Advanced as part of a long-term, consistent skin protocol, not as a one-time or medical intervention.</p>
--	---

Method overview

Parameter	Description
Participant	Female participant, Skin Type 2 / fair skin. Personal identifiers are excluded from this public report.
Assessment system	Aura Face Analysis - skin scoring across wrinkles, texture, brown spots, red areas and pores.
Dates	Baseline: 13.10.2025. Follow-up: 20.02.2026.
Product context	LumiTech+ Advanced LED Face Mask used as part of a regular home skincare protocol.
Analysis type	Single-subject observational follow-up. Not a controlled clinical trial.

<p>Regulatory language note</p> <p>Internal Skin Follow-Up, Observational Skin Analysis Report.</p>
--

Observed results

All values are interpreted in the context of the Aura 1-5 scoring system, where a lower value indicates a more favorable result for the respective parameter. The relative skin score uses +/- comparison against a reference group.

Skin metric	Baseline	Follow-up	Change	Interpretation
Texture	1.6	1.3	-0.3	Directional improvement
Brown spots	2.7	2.3	-0.4	Directional improvement
Wrinkles	1.2	1.1	-0.1	Small directional improvement
Red areas	1.2	1.2	0.0	Stable
Pores	2.1	2.1	0.0	Stable
Relative skin score	0.0	-0.2	-0.2	Improved vs reference group

Direction of change on 1-5 scale

Metric	1-5 score movement	Change
Texture	1.6 -> 1.3	-0.3
Brown spots	2.7 -> 2.3	-0.4
Wrinkles	1.2 -> 1.1	-0.1
Red areas	1.2 -> 1.2	0.0
Pores	2.1 -> 2.1	0.0

Main interpretation


The most notable observed change was seen in brown spots and texture. Stable red areas and pores are also relevant, as no worsening was observed within the follow-up period.

These results are observational and should be viewed as part of an overall skin routine, not as isolated proof of the LED mask effect alone.

Visual endpoint excerpts


The excerpts below are from the original Aura Face Analysis pages for the respective parameter. They show regional and overall scoring maps for baseline and follow-up.

Texture - skin surface evenness

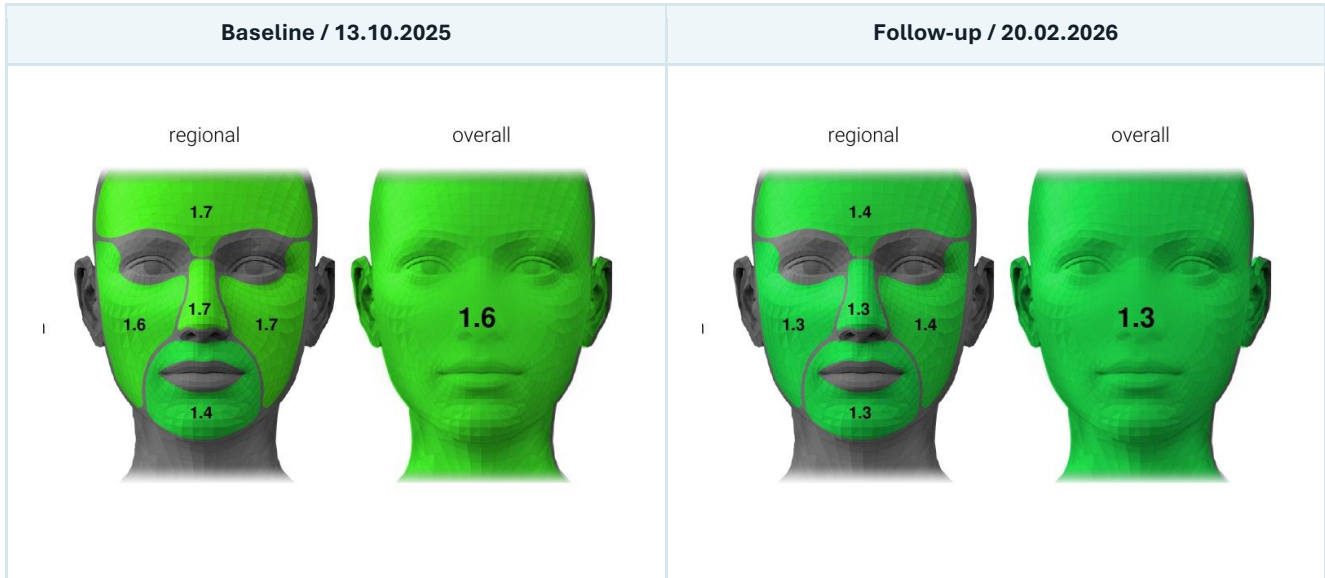Texture score decreased from 1.6 to 1.3, indicating a more favorable measured skin texture score on the Aura 1-5 scale.

Brown spots - pigmentation-related score

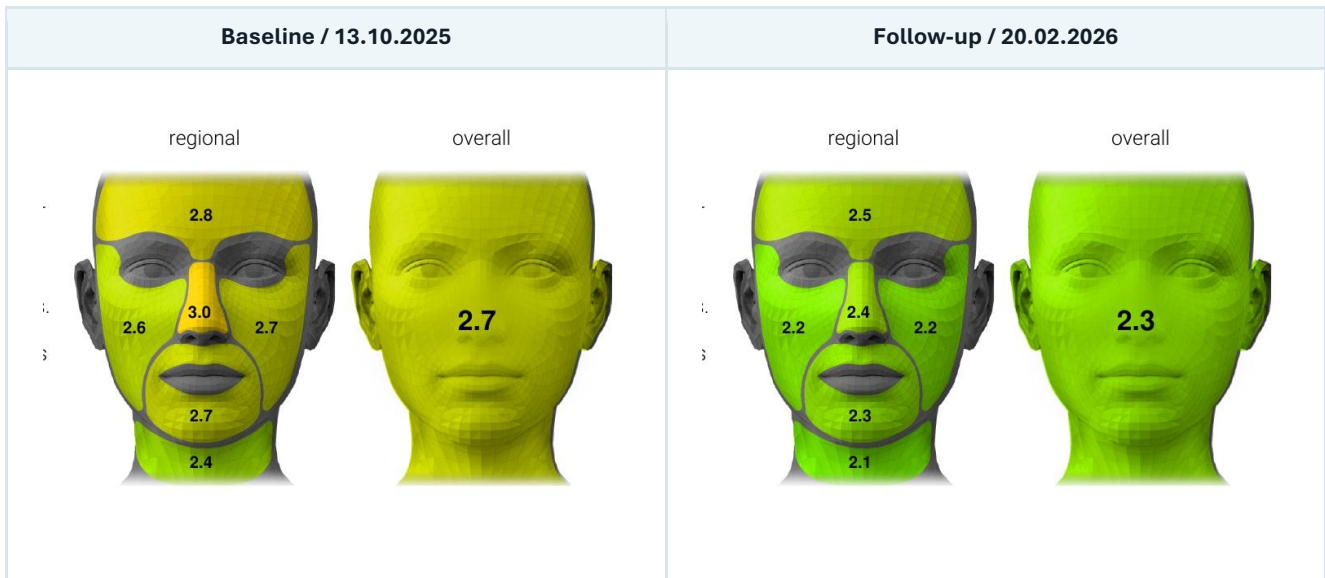
Brown spots score decreased from 2.7 to 2.3. This is a positive directional change, but pigmentation can also be influenced by sun exposure, seasonality, sunscreen consistency and topical skincare actives.

Secondary endpoint excerpts


These parameters show a small positive change or stability.

Wrinkles

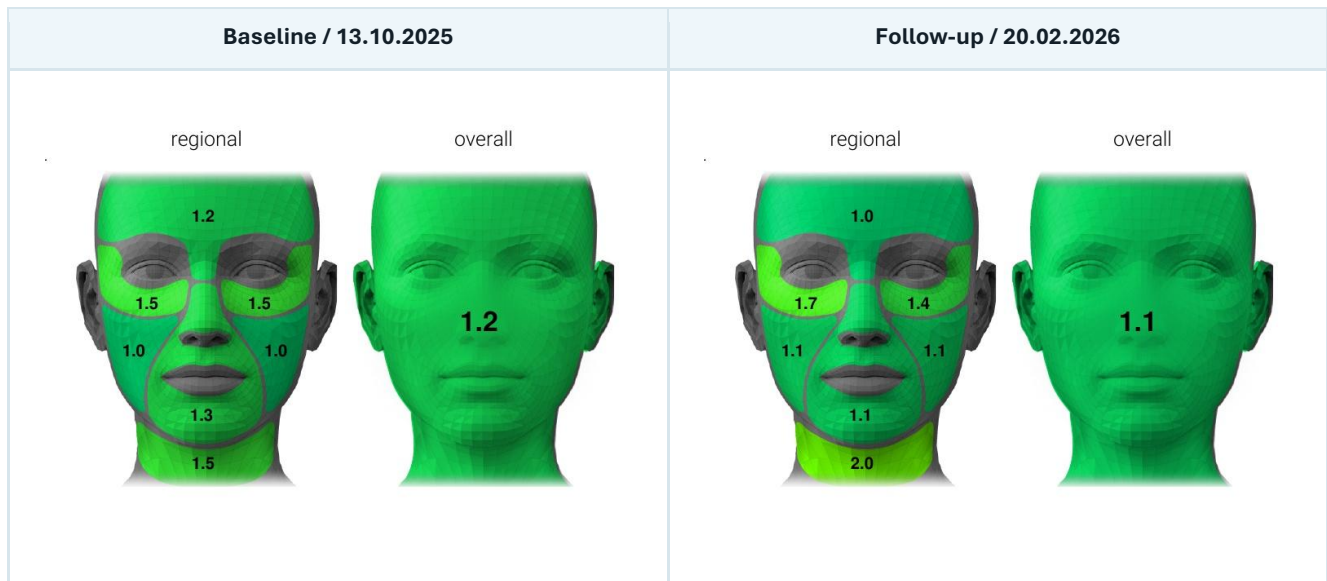Wrinkles score decreased from 1.2 to 1.1 - small directional improvement.

Red areas

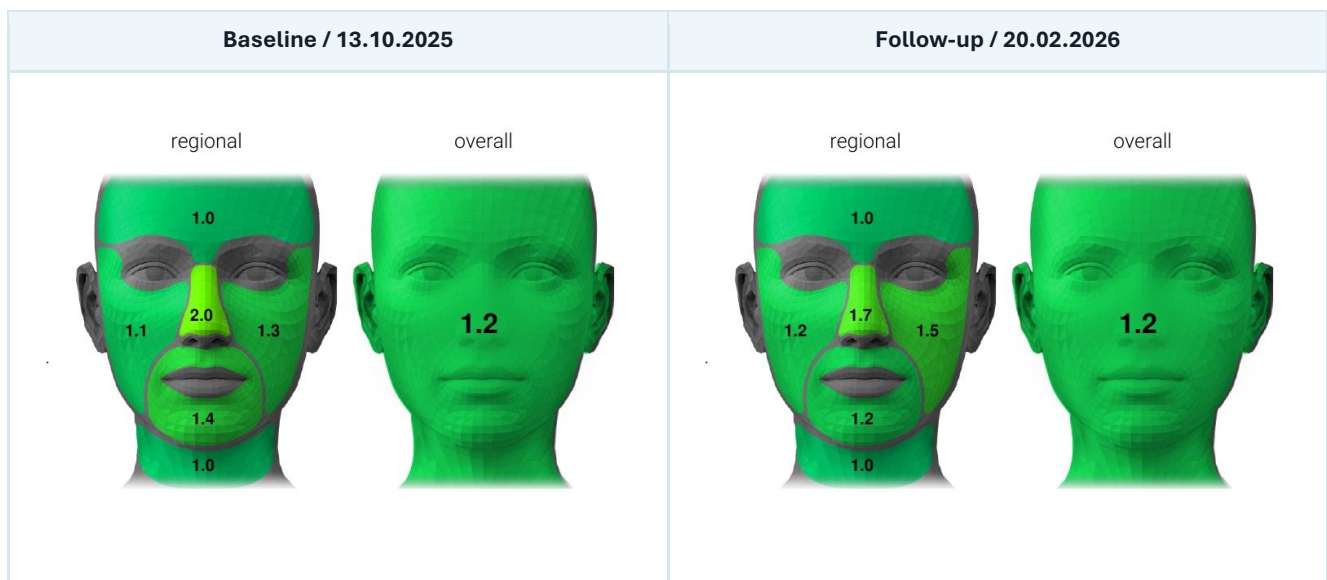
Red areas remained stable at 1.2. Stability is relevant in a home-use device context.

Pores - stability endpoint

Pores remained stable within the observed period. This is a moderate but relevant observation, because no worsening of the parameter was measured.

Pores

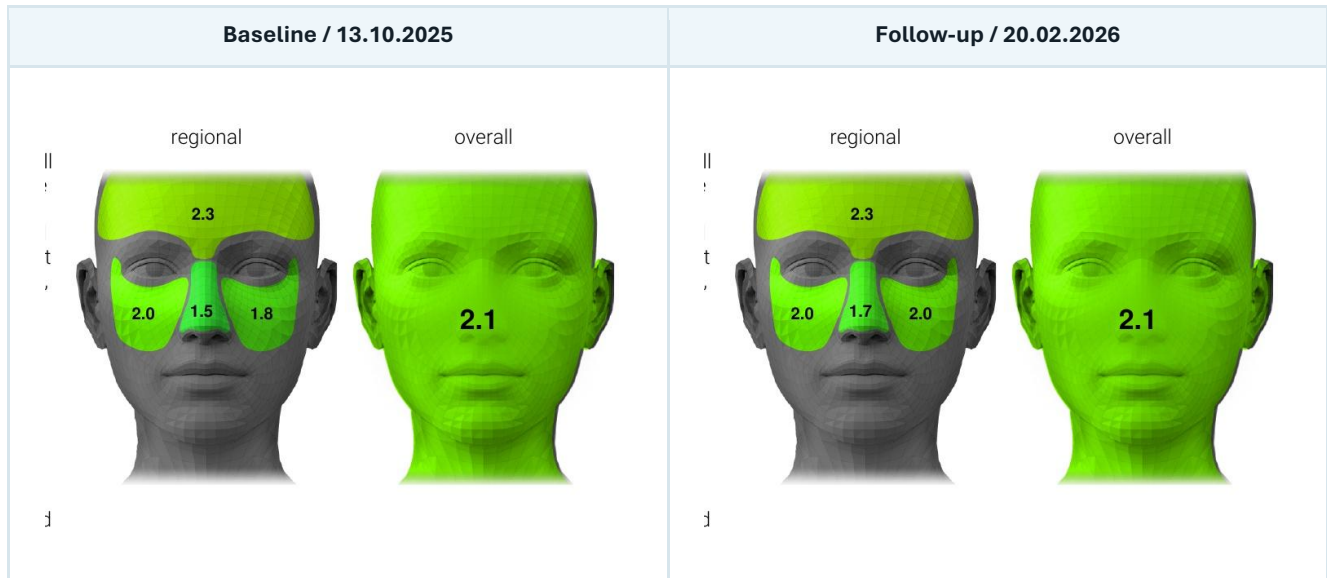
Pores remained stable at 2.1. No worsening was observed in this follow-up.

Stability interpretation

Stable results for red areas and pores.

Aura scoring overview

The Aura report provides two complementary scoring views: absolute 1-5 scoring and relative +/- scoring compared to a reference group with the same age, gender and skin type. The excerpts below are included to show the original scoring context used for this follow-up.

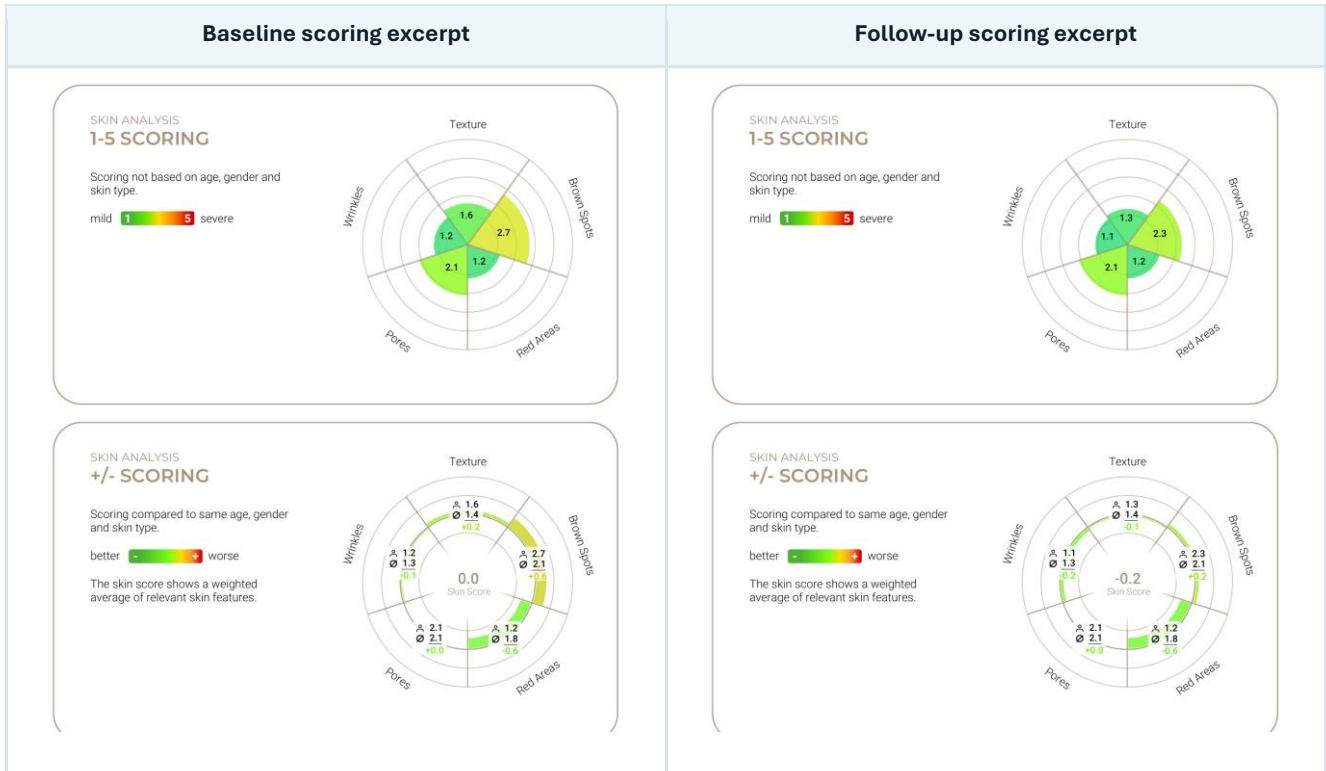Source: Aura Face Analysis reports, capture dates 13.10.2025 and 20.02.2026. The original documents include personal identifiers; those are intentionally excluded from this public-facing AL+ version.

Technology rationale

LumiTech+ Advanced is a home-use LED face mask positioned as a cosmetic/wellness device. It is intended for regular, non-invasive use as part of a skin routine and does not replace medical consultation or treatment.

Wavelength zone	Cosmetic rationale
Red light around 630 nm	Supports cosmetic care focused on skin vitality, comfort and visible radiance.
Yellow light around 590 nm	Used in cosmetic LED routines focused on tone comfort and a calmer-looking complexion.
Blue light around 460 nm	Used in routines targeting the appearance of oilier, blemish-prone skin.
Near-infrared around 850 nm	Used in LED routines focused on deeper cosmetic support and long-term skin appearance.

Positioning statement

The observed results are consistent with the intended cosmetic positioning of LumiTech+ Advanced: long-term support for a more even, balanced and cared-for skin appearance when used regularly and responsibly.

Limitations and responsible interpretation

- Single-subject observational follow-up.
- No placebo or control group.
- No randomization or blinding.
- The LED mask was part of a broader skincare routine.
- Seasonality, sun exposure, sunscreen habits, sleep, hormones and topical products may influence the measured skin parameters.
- Results are individual and should not be presented as guaranteed or universal.

Public claim boundary

Internal observational follow-up reported improvement in selected Aura skin analysis parameters.

Conclusion

Between 13.10.2025 and 20.02.2026, Aura Face Analysis reported positive directional change in texture score, brown spots score, wrinkles score and relative skin score, with stable red areas and pores. These data support the role of LumiTech+ Advanced as a premium LED care device within a consistent, long-term skin routine.

This document is intended for transparency, product education and internal evidence-building. It is not medical advice and does not diagnose, prevent, treat or cure disease.

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AL+ Cosmetics / ALMA DERM Ltd. / alpluscosmetics.com